INCREIBLE<sup>®</sup>

DEVICES

# ASTRODOME<sup>™</sup>

G L A S S F A C I A L


## OxygenCeuticals Korean Glass Facial Astrodome Protocol

Designed for superior facial treatments, the **OxygenCeuticals Korean Glass Astrodome Facial<sup>™</sup>** Machine offers a range of skin-enhancing solutions. This Protocol Guide offers post-care and aesthetic solutions for Medical Aesthetic Professionals.

 HEALTH CANADA LICENCED

# Astrodome Facial Postcare System
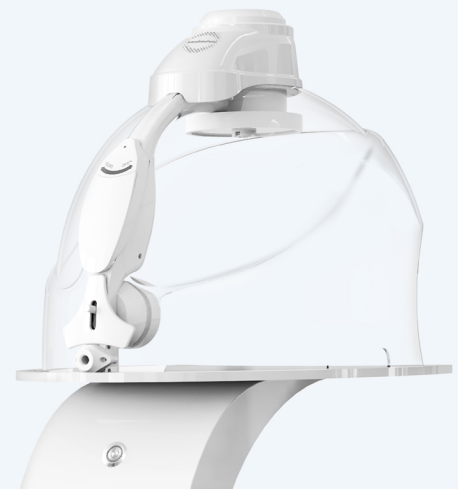
## 01 Face & Scalp Treatment

- Dual LED Panels to treat face and scalp at the same time
- Paneles duales para tto de cara y cuero cabelludo al mismo tiempo


## 02 Oxygen & Anion Therapy

- Oxygen & Anion mist infusion with targeted OxygenCeuticals solutions
- Spray infusion de anion y oxigeno con objetivo OxygenCeuticals soluciones

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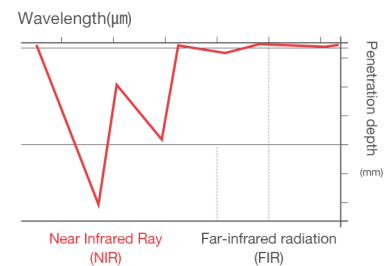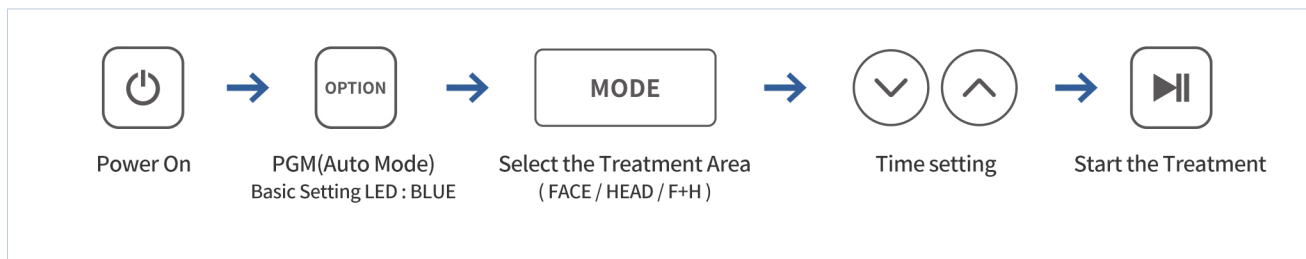
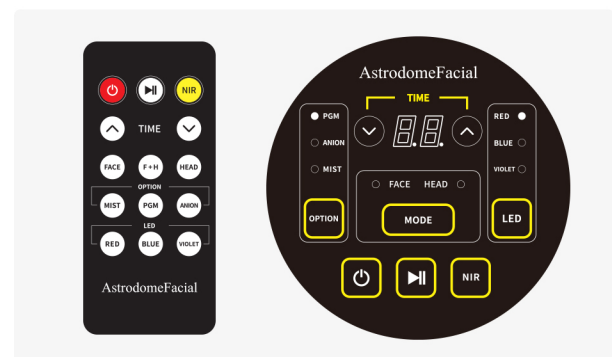
## 03 Oxygen Infusion

- Infuse nutrients based on skin types
  - preparado de nutrientes basados en tipos de piel
- [ Use ] OxygenCeuticals Formula 5~10ml

## 04 Oxygen Mask ( I )

- Apply HA Gel+Toning Gel
- Aplicar Gel y Tónico HA Gel

## 05 Inject the Cryogenic Activator

- Inject 6 to 8cc Cryogenic Activator
- Inyectar 6 a 8cc Cryogenic Activator

## 06 AstrodomeFacial Treatment

- Select suitable program
- Seleccionar Programa ajustable

## 07 Oxygen Mask ( II )

- CryoGel Velvet Mask **or** Oxygen Mask(RT/CT/VC)
- Mascara criogel terciopelo u mascara oxigeno (RT/CT/VC)

## 08 Scalp Scaling + Finishing


- Select suitable finishing products
  - Selección final adecuada productos
- [ Use ] Cryogenic Activator

# AstrodomeFacial Automatic Treatment Program

## Acne

BLUE

- Relieve skin troubles ( acne, excessive sebum )
- Alivia problemas de piel ( acne, sebo excesivo )

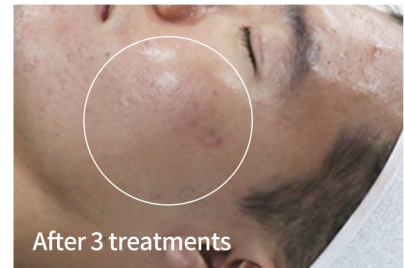**Finish :** Toning Gel → Herbal Fluid → Purifying Wintergreen Serum  
→ Purifying Wintergreen Cream or Cell Couperose Cream

## Brightening

RED

- Brightening and firming
- Aclaramiento y firmeza


**Finish :** Toning Gel → Fermented Fluid → Moisture Aqua Serum → PP Cream or Hydrating Cream

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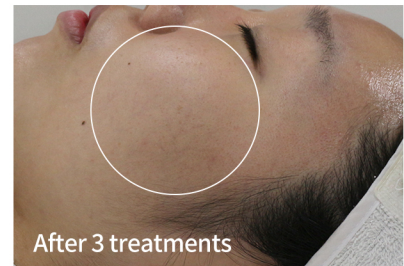 [incrediblelaser.com](http://incrediblelaser.com)

 [hello@incrediblelaser.com](mailto:hello@incrediblelaser.com)

 [facebook.com/incrediblelaser](https://facebook.com/incrediblelaser)

 [instagram.com/incrediblelaser](https://instagram.com/incrediblelaser)

 [tiktok.com/@incredibledevices](https://tiktok.com/@incredibledevices)

 [dailymotion.com/incredibledevices](https://dailymotion.com/incredibledevices)

 **HEALTH CANADA LICENCED Medical Devices**

**All products are compliant with Health Canada Safety regulations.**